



U.S. Patents 6,571,881 &
7,533,733

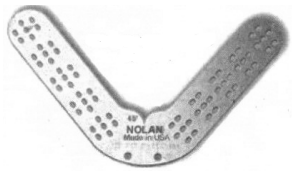

NOLAN HOOF PLATE COMPANY

Solutions for Hoof Problems

We want to help!

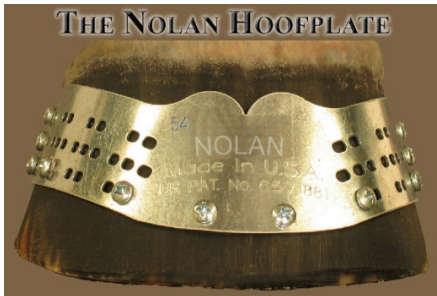

*The only patented product
shown to assist
therapy of a
variety of
hoof problems.*



**HOOF CRACKS, FLARES, THIN SOLES, POOR WALL STRUCTURE,
LAMINITIS AND MORE. FARRIER CERTIFICATION AVAILABLE.**

NOLAN HOOF PLATE COMPANY
2417 Jericho Turnpike, #333
New Hyde Park, NY 11040

INFO: (614) 565-2810
PHONE ORDERS: (888) 600-0360
FAX ORDERS: (888) 600-0310

**OVER 100 VETS AND FARRIERS CERTIFIED AND APPLYING
THE NOLAN HOOF PLATE ACROSS NORTH AMERICA**

NOLAN HOOF PLATES

Best uses: Assists horses in recovery from many hoof problems to including wall cracks, thin and flat soles, laminitis, and much more.

How it works: The Nolan Hoof Plate and horseshoe work together as an integral unit to provide external support to the hoof wall from the toe to both caudal quarters. This unit then functions to create hoof stability, with greatly increased wall strength. This stronger wall no longer over-flexes which allows for greater blood flow throughout the vessels of the hoof capsule. The improved circulation then stimulates the laminae to form horn at an accelerated rate.

How to use it: Correctly trim and prepare the hoof for shoe and plate application. Since the plate and shoe will remain in place for up to 16 weeks, preparation is crucial. The hoof angle is measured and the proper plate is selected, formed and applied with the shoe. The plate is then covered with an appropriate adhesive to further enhance, stabilize and protect the application. Each Nolan Hoof Plate Kit contains complete installation instructions.

Comments: Professional farrier Mike Nolan of Columbus, Ohio, originally created the Nolan Hoof Plate to help keep horseshoes on working draft horses, particularly in difficult terrain where rocks and mud often cause shoes to be pulled off. Current on-going clinical research shows additional benefits of the Nolan Hoof Plate when applied to assist a variety of hoof maladies facing farriers today, providing greater hoof wall stability and improved blood flow. The Nolan Hoof Plate has also been used as a cast to provide excellent stability to facilitate healing of coffin bone fractures.